Upgrade · The Home Creator Way

LAB DISTILLATION APPARATUS

USER MANUAL

MODEL:YHEO-3

MODEL:YHEO-3

This is the original instruction, please read all manual instructions carefully before operating. VEVOR reserves a clear interpretation of our user manual. The appearance of the product shall be subject to the product you received. Please forgive us that we won't inform you again if there are any technology or software updates on our product.

SAFETY INSTRUCTIONS

WARNING:

Read this material before using this product. Failure to do so can

result in serious injury.

Assembly precautions

1. Assemble only according to these instructions. Improper assembly can create hazards.
2. Wear safety goggles and duty work gloves during assembly.
3. Keep the assembly area clean and well-lit.
4. Keep bystanders out of the area during assembly.
5. Do not assemble when tired or when under the influence of alcohol, drugs or medication.
6. The product capabilities apply to properly and completely assembled products only.
7. Assemble on a flat, level, hard and smooth surface capable of safely supporting the Lab Distillation Apparatus.
8. For additional information regarding the parts listed in the following pages, please refer to the Assembly Diagram of this manual. Unwrap and separate all parts in a clean work area.

Use precautions

1. DO NOT SIT OR STAND ON THIS ITEM.
2. This product is not a toy. Do not allow children to play with or near this item.
3. Do not exceed weight capacities.
4. Use only on a flat, level, hard and smooth surface capable of safely supporting a fully loaded Lab Distillation Apparatus.
5. Use as intended only.
6. Inspect before every use; do not use if parts are loose or damaged.
7. The product complies with relevant standards and regulations.
Installation of this device must comply with all applicable regulations, and it should only be operated in well-ventilated areas. Please read the user manual before installing and using this equipment.
8. This appliance is not intended for use by persons (including children) with reduced physical, sensory or mental capabilities, or lack of experience and knowledge, unless they have been given supervision or

instruction concerning use of the appliance by a person responsible for their safety.

9. Children shall not play with the appliance
10. This appliance shall not be used by children. Keep the appliance and its cord out of reach of children.
11. This appliance is intended to be used in household and similar applications such as:
 - Offices environments;
 - farm houses;
12. Please properly keep product. Check whether the product is abnormal before use. If it is abnormal, please stop using it
13. Please use the product according to the intended use, warn of potential injury from misuse.

SAVE THIS MANUAL

INSTRUCTIONS


Thank you very much for choosing this model. Glass fittings should always be handled with care. Please read all of the instructions before using it. The information will help you achieve the best possible results.

Parts name	Number of parts	Description
Flask	1	500ml
Coil condenser	1	/
Pear-shaped Separatory funnel	1	125ml
Conical flask	1	100ml
Glass stopper	2	19#
beaker	1	100ml
Glue head dropper	1	Red glue head
Alcohol burner	1	150ml
thermometer	1	100℃
Rubber plug	4	/

Latex tubing	1	2m
Vacuum grease	1	Vaseline
PH test paper	1	/
Flask brush	1	Small size
Iron support stand	2	With 40cm supporting pole
Flask clamp	2	/
Cross clamp	2	/
Horn tube	1	/
Circular ring	2	1 of 5.1cm, 1 of 7.3cm
Terracotta net	1	/

OPERATION


1. Assemble the iron support stand;

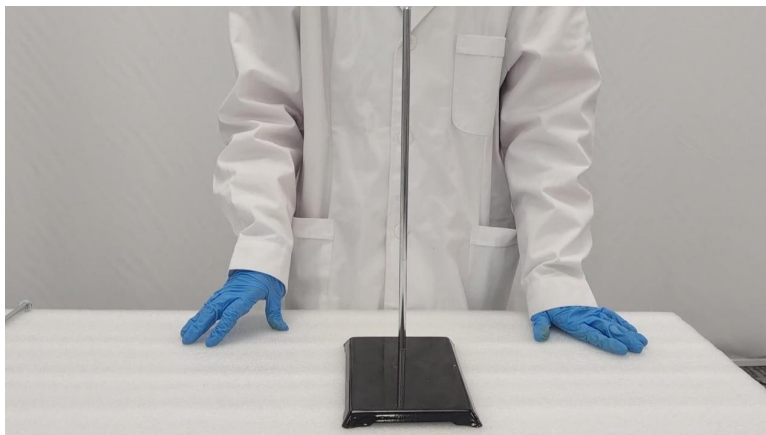
2. Place alcohol burner;


3. Install the circular ring;

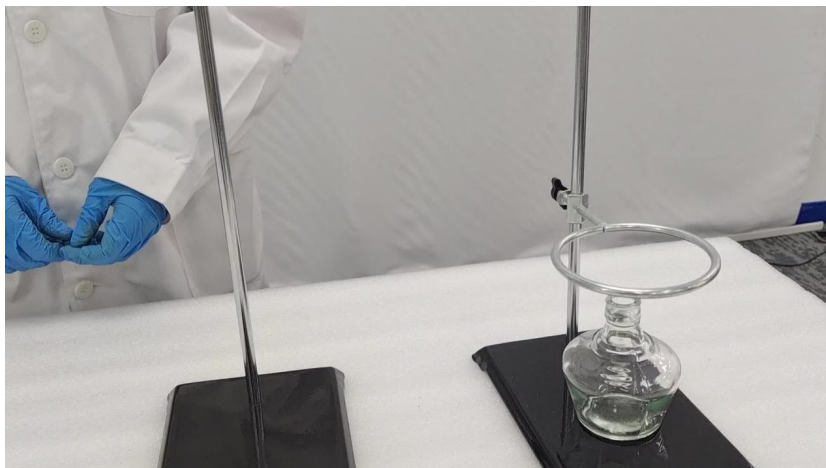
4. Place the terracotta net on the circular ring , and install the cross clamp and flask clamp at the appropriate height;

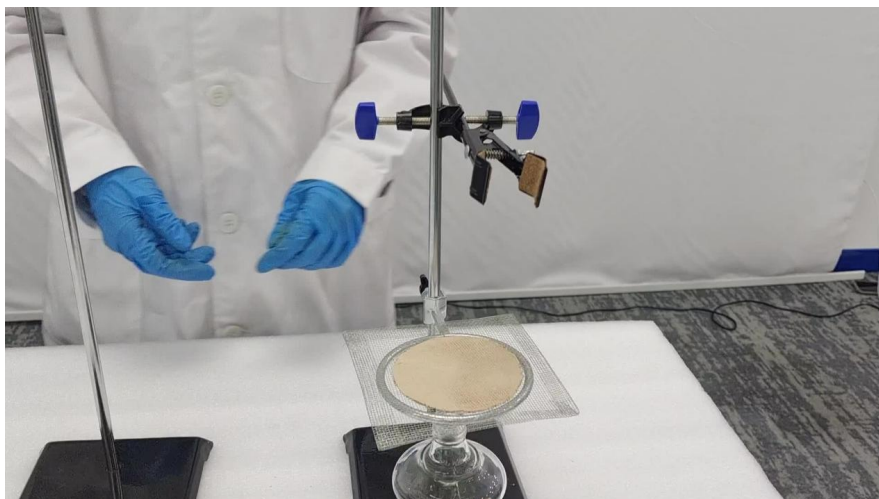
5. Install the flask, and install the cross clamp at the appropriate height of the other iron support stand;

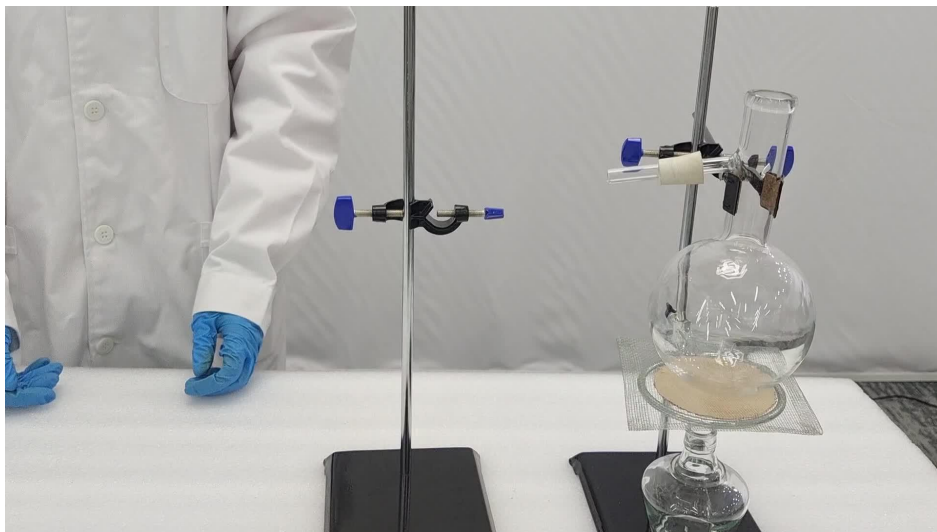
6. Install the coil condenser and fix it with the flask clamp;

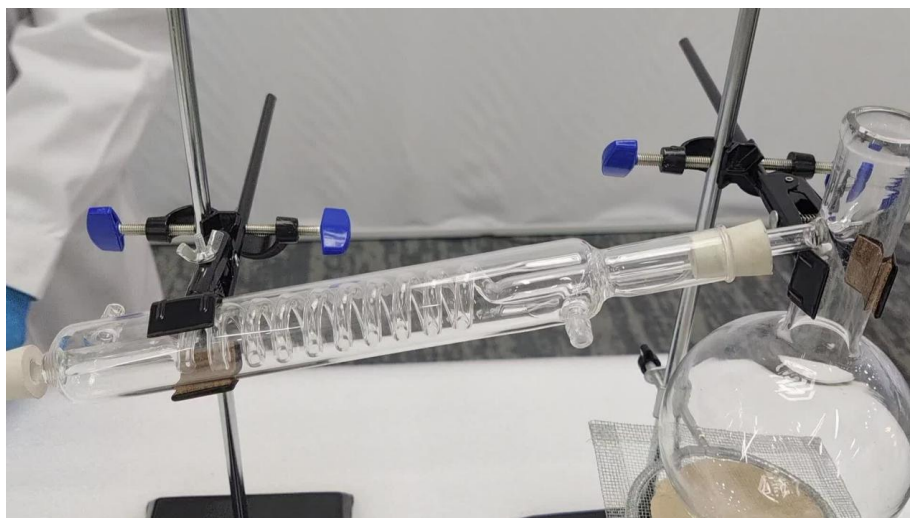
7. Install the horn tube and place the conical flask under the outlet of horn tube;

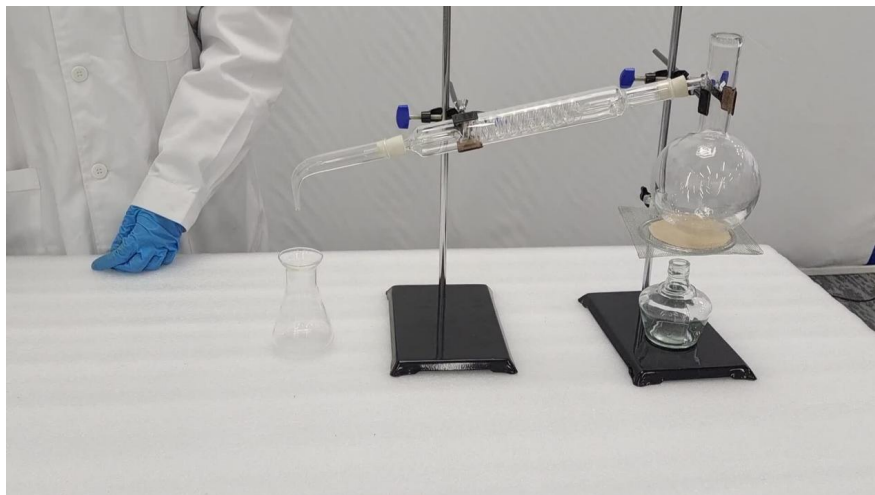8. Install the circular ring at the appropriate height of the iron support stand, and place the pear-shaped separatory funnel on it, then place the beaker below the outlet of the pear-shaped separatory funnel;

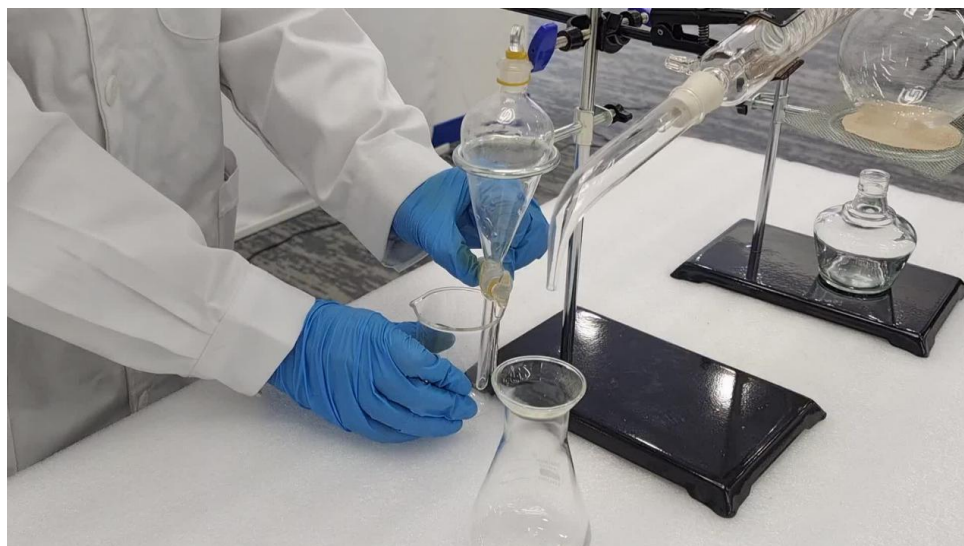9. Screw the thermometer into the rubber plug and install the rubber plug in the mouth of the flask.

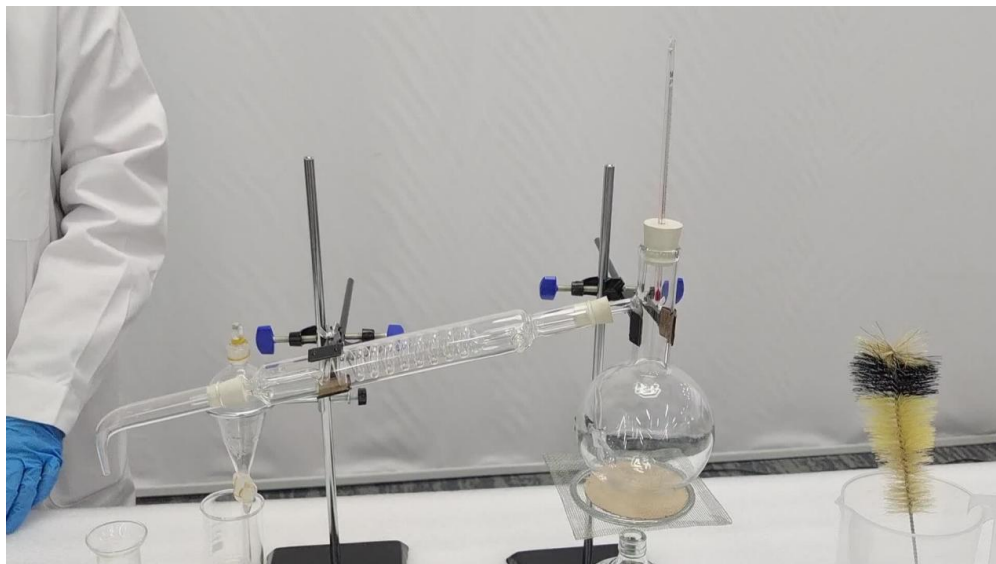PRODUCT PARAMETER

Model	YHEO-3
Heating Method	Alcohol Burner
Flask Type	Distilling Flask

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